

Intelligent 5-Axis Dry Dental Milling Machine

BSM-500D/BSM-500DS Compressor-Free Version

Key Features



Precise

- Digital servo system with high resolution
- Guaranteed milling accuracy: ±0.005mm
- Auto-calibration ensures consistent precision



Efficient

- 60,000 RPM German spindle for high-speed milling
- 1/3 holder boosts output by 20% per block
- Milling speed:
 - Coping: 11 mins(Zirconia), 5 mins(Wax)
 - Crown: 16 mins(Zirconia), 6 mins(Wax)
 - Splint: 45 mins(PMMA/PEEK)

Above data were obtained based on the BSM-500D milling of Aconia HT+ 14mm



Reliable

- Industrial-grade aerometal gantry structure
- Advanced dust-proof system protects critical components
- Integrated ionizer neutralizes static for effortless cleaning



Intelligent

- Automatic tool changer
- Power-off protection & break-point resume
- Easy operation with intuitive LCD touch screen
- Automatic working-suspension function
- Intelligent working status indication light



BSM-500D, with a strong & reliable gantry structure design and 5-axis simultaneous linkage, features excellent reliability, intelligent processing, and efficient, precise milling at 5µm milling accuracy. With an automatic working-suspension function in case of accidental door opening during processing, it ensures operator safety to avoid physical injury. The built-in working status indicator provides intelligent monitoring of the milling operation. Equipped with a compact yet powerful spindle, it is an excellent choice for both dental labs and chairside practices.

Millable Materials



Zirconia



PMMA



Wax



Wood



PEEK

Millable Restorations



Full contour crown



Coping



Inlay/onlay



Veneer



Crown bridge



Zirconia abutment



Zirconia implant bridge



Non-metal frame

Technical Parameters

Width*depth*height	700mm*480mm*645mm	Weight	90kg		
Number of axis	5	Motor	AC servo		
Spindle power	800W	Spindle speed	60000RPM		
Built-in-burs	6	Rotation angle	A axis: 360°, B axis: -30°~ + 30°		
Tool Type	Zirconia End Mill	φ4mm*1mm			
	Zirconia Ball Mill	φ4mm*2mm	φ4mm*1mm	φ4mm*0.6mm	
	Resin Mill	φ4mm*2mm	φ4mm*1mm		
Voltage/Power	220V/1KW	Fields of application	Dry milling		
Milling accuracy	±0.005mm	External Control Linkage	Vaccum		
Safety function	Automatic working-suspension function	CAM Software	Open system		
Remote control	√	Monitor	8"LCD touch screen		

BSM-500DS

Compressor-Free Version

Built for Compact Spaces & Quiet Operation



Efficiency Upgrade

- **Lower maintenance:** No compressor upkeep required
- **Energy optimized:** Direct-drive cuts production energy use



Space & Noise Control

- **Compact design saves space:** Ideal for clinics/small workshops
- **Ultra-low noise:** Vibration-free operation reduces sound output



Easy Operation

- **Plug-and-play setup:** No complex air lines needed
- **Instant production:** Start producing instantly at power-on

Feature Comparison

BSM-500DS



Pneumatic-free design

0 compressor upkeep

Clinic-quiet & desktop-ready

BSM-500D

Compressor + air lines

Annual maintenance cost

Standard noise output

CAD/CAM Digital Processing



About Besmile

Besmile is a global digital solutions provider for restorative and implant dentistry.

We offer integrated systems—premium CAD/CAM materials, advanced equipment, and precision implant systems—designed to streamline workflows and ensure reliable results.

All core products are developed and manufactured in-house, ensuring consistent quality and continuous innovation.

Trusted by over 1,000 partners in more than 100 countries, we empower dental professionals to create confident, lasting smiles.

Technology creates the best smile.



Chengdu Besmile Biotechnology Co.,Ltd.

www.bsmdental.com

T: +86-28-85317108

E: info@cdbesmile.com


Besmile Dental America Inc.

20311 Valley Blvd, Suite#1, Walnut, CA91789

T: (626)921-5798

E: sales@bsmdentalus.com

Contact Us

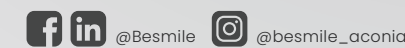

Intelligent 5-Axis Dental Milling Machine

BSM-500D BSM-500DS

